



## VG55 VENTILATOR



**Exclusive Distributor Aeonmed & Heyer For Sulawesi**

**PT KARYA PRATAMA**


Daan Mogot Arcadia Blok F1 No. 8A-8B Jl. Daan mogot Km 21,  
Kec. Batuaceper, Kel. Batuaceper, Kota Tangerang, Banten 15122  
Telepon : (021) 5595 8085, Fax : (021) 5595 8125




# VG55 Ventilator

## Start Your Natural Breath with VG55


VG55 is specially designed for Non-Invasive Ventilation (NIV). Research has shown that early application of NIV can effectively decrease intubation rates and hospitalization length. The application of NIV can also prevent ventilator associated pneumonia (VAP), as well as promote respiratory care. VG55 adopts advanced turbine technology, provides NIV with ICU-grade monitoring.

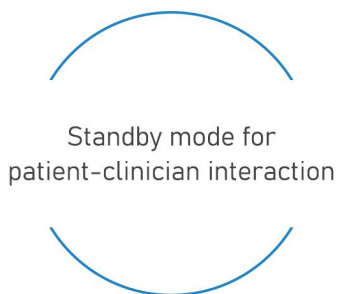

Standby mode for patient-clinician interaction

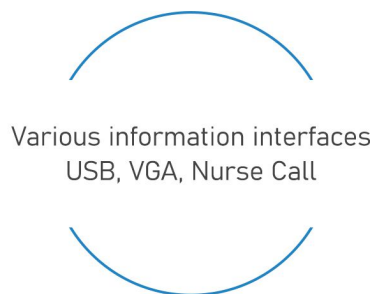

Various information interfaces  
USB, VGA, Nurse Call

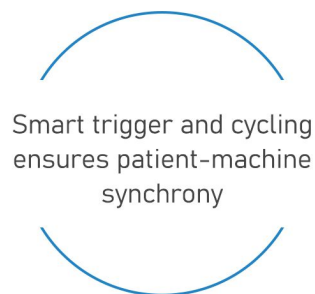

Smart trigger and cycling ensures patient-machine synchrony



Alarm reset helps maximize patient safety while reducing excessive alarms



Exclusive Distributor Aeonmed & Heyer For Sulawesi

**PT KARYA PRATAMA**

Daan Mogot Arcadia Blok F1 No. 8A-8B Jl. Daan mogot Km 21,  
Kec. Batuaceper, Kel. Batuaceper, Kota Tangerang, Banten 15122  
Telepon : (021) 5595 8085, Fax : (021) 5595 8125

**KONTAK DIVISI E-CAT :**

0812 9351 3876 (Kamal)  
0822 1148 1678 (Endah)



**15** Seconds

## Start Your Natural Breath within 15 Seconds

- Auto-calibration is performed within 15 seconds, safely guaranteeing the ventilator's performance.
- Quick setup with guide-way system and intuitive interface.
- Smart trigger combines pressure and flow triggers, highly improving the patient-machine synchrony, following the patient to breathe in a more natural manner

**5** Minutes

## Become a Respiratory Therapist in 5 Minutes

- With 15,6" intelligent touchscreen, the VG55 displays complete ventilation monitoring in a single screen.
- Simplifies the patient management with user-friendly interface and guide-way information.
- 3000 pieces event/alarm logs provide clear information to ensure continuity of patient care.

**in** Seconds

## Get What You Need in Seconds

- Two minutes of emergency oxygen therapy provides 100% oxygen for whenever and wherever needed.
- Extended SpO<sub>2</sub>, PR module fulfill clinical requirements.
- High-Flow Therapy up to 60L/min.



**Exclusive Distributor Aeonmed & Heyer For Sulawesi**

**PT KARYA PRATAMA**

Daan Mogot Arcadia Blok F1 No. 8A-8B Jl. Daan mogot Km 21,  
Kec. Batuaceper, Kel. Batuaceper, Kota Tangerang, Banten 15122  
Telepon : (021) 5595 8085, Fax : (021) 5595 8125


**KONTAK DIVISI E-CAT :**

0812 9351 3876 (Kamal)  
0822 1148 1678 (Endah)

# VG55 Ventilator

## High-Flow Therapy

- An air/oxygen blender, allowing for 21% to 100% FiO<sub>2</sub>, generates up to 60L/min flow.

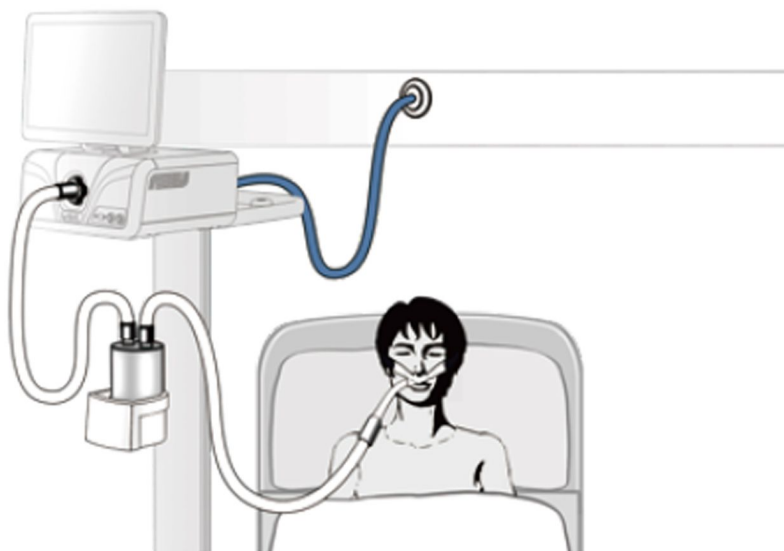

## Complete Ventilation Modes

- Volume Guaranteed Pressure Support Ventilation (VGPS) assures tidal volume with auto adjusting pressure support, decreasing the workload for clinicians while ensuring patient safety
- CPAP with Crelief option offers two levels of flow-based expiratory pressure relief, which improves patient comfort and treatment acceptance.



Exclusive Distributor Aeonmed & Heyer For Sulawesi

**PT KARYA PRATAMA**

Daan Mogot Arcadia Blok F1 No. 8A-8B Jl. Daan mogot Km 21,  
Kec. Batuaceper, Kel. Batuaceper, Kota Tangerang, Banten 15122  
Telepon : (021) 5595 8085, Fax : (021) 5595 8125

**KONTAK DIVISI E-CAT :**

0812 9351 3876 (Kamal)  
0822 1148 1678 (Endah)